




## INTEGRATED PLATFORM FOR DRUG DESIGN AND DISCOVERY

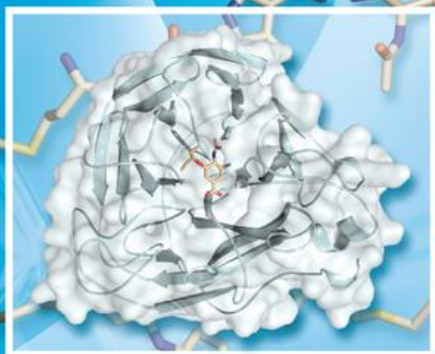

### CUSTOMIZED WORKFLOWS FOR DRUG DISCOVERY

Automated docking using *FITTED*

*Filter by descriptors*

*Convert 2D to 3D*

*Prepare protein - pdb to mol2*

*Extract representative library*

*Create combinatorial library*

*Search for analogues*

*Easy integration of 3rd party programs*

[www.molecularforecaster.com](http://www.molecularforecaster.com)

Molecular Forecaster Inc.

FORECASTER Platform 1.2 – User Guide

13 June 2014

## TABLE OF CONTENTS

Introduction .....	5
The Regular User Interface.....	6
General Comments – Proper Usage.....	6
Minimum requirements for clients.....	6
Login Interface.....	6
Home Interface .....	6
Working with Forecaster.....	7
The Sections .....	7
Workflow Manager .....	7
Adding a new workflow.....	8
Building a workflow diagram .....	8
Jobs Manager .....	9
Adding a new job.....	9
Editing the parameters and selecting the files to process.....	10
Executing and viewing jobs .....	13
File preview .....	13
File Manager.....	14
Folder - My Files .....	14
Folder - My Proteins.....	14
Job-Specific Folders.....	15
Files and actions .....	15
Sketcher.....	16
Saving and Loading files .....	16
Reactions Manager .....	17
Creating reactions .....	17

Administrator Interface - Settings.....	19
User Manager .....	19
User roles and access .....	19
Adding new users .....	19
Plugins Manager .....	20
Adding new plugins .....	20
Adding parameters.....	22
Adding new plugins from XML .....	23
Configurations Manager .....	24
Default settings .....	24
Save Database to Fixtures .....	24
Actions.....	25
Structure-based.....	25
Dock ligand(s) using Fitted .....	25
Predict som using impacts .....	26
Proteins/Nucleic acids.....	27
Make proteins similar.....	27
Superpose protein structures .....	28
Prepare protein – pdb to mol2 .....	28
Setup protein for docking .....	30
Ligand-based .....	30
Add descriptors .....	30
Filter by descriptors.....	31
Combinatorial Library from 2D structures .....	31
Search for analogues.....	32
Extract representative library .....	32

Combinatorial Library from 3D structures .....	33
Add fragments to molecules .....	34
Pharmacophore identification and search.....	34
Substructure search .....	35
Small molecules.....	35
Setup ligand(s) for docking.....	35
Add hydrogen to 2D structures.....	35
Convert 2d to 3d.....	35
Clean structure geometry .....	36
Setup reactants for combichem in 3D .....	36
Structure optimization .....	36
Other .....	36
Convert file formats .....	36
Function linker.....	36
References.....	37

## INTRODUCTION

The programs included in the FORECASTER Platform have been developed within the Moitessier group at McGill University.

Team leader: Prof. Nicolas Moitessier

Development and support team leader: Dr. Eric Therrien

<http://www.molecularforecaster.com/support.html>

**Please refer to the FORECASTER platform by citing the following publication:**

Therrien E., Englebienne P., Arrowsmith A.G., Mendoza-Sanchez R., Corbeil C.R., Weill N., Campagna-Slater V., Moitessier N. Integrating medicinal chemistry, organic/combinatorial chemistry, and computational chemistry for the discovery of selective estrogen receptor modulators with FORECASTER, a novel platform for drug discovery "Journal of Chemical Information and Modeling" **2012**, 52, 1, 210-224

**Please refer to the FITTED program by citing the following publication:**

Corbeil C.R., Englebienne P., Moitessier N. Docking ligands into flexible and solvated macromolecules. 1. Development and validation of FITTED 1.0 "Journal of Chemical Information and Modeling" **2007**, 47, 2, 435-449

**Please refer to the IMPACTS program by citing the following publication:**

Campagna-Slater V., Pottel J., Therrien E., Cantin L.-D., Moitessier N. Development of a computational tool to rival experts in the prediction of sites of metabolism of xenobiotics by P450s "Journal of Chemical Information and Modeling" **2012**, 52, 9, 2471-2483.

## THE REGULAR USER INTERFACE

### GENERAL COMMENTS – PROPER USAGE

When using the platform, it is important to **avoid any white space in filenames**. The platform might recognize the file, but the programs will fail to read the correct filename. Use “\_” or “-” instead.

### MINIMUM REQUIREMENTS FOR CLIENTS

FORECASTER has been tested on Internet Explorer 7 and higher, Firefox 3.6 and higher, Chrome, and Safari 4 and higher. Other browsers including Internet Explorer 6 are not fully compatible with the interface FORECASTER. However, they can still be used but some logos, fonts and other visual details may be affected.

### LOGIN INTERFACE

When loaded in a browser, FORECASTER will display the login interface. User inputs login name and password to log into the system. The default account is “admin” with password “fitted” (without the quotes).

LOGIN	<input type="text"/>	<div>LOGIN</div>
PASSWORD	<input type="password"/>	


### HOME INTERFACE

When logged in, the user will see the home page. This page presents the main features of the application.

The screenshot displays the FORECASTER Molecular web interface. At the top, the logo "MOLECULAR FORECASTER" is visible alongside the tagline "Discovering Molecules". A "Guest User" button is in the top right corner. Below the header is a navigation bar with links: Home, Contact Us, FAQ, and LOGOUT. A secondary navigation bar contains links to various modules: HOME (highlighted), WORKFLOWS MANAGER, JOBS MANAGER, FILE MANAGER, SKETCHER, REACTIONS, and SETTINGS. The main content area is titled "Forecaster - Drug Discovery Suite" and features six interactive tiles, each with a molecular icon and a description of a function: "Dock ligand(s) using FITTED", "Setup protein for docking", "Filter by descriptors", "Add descriptors", "Extract representative library", and "Create combinatorial library".


### THE SECTIONS

When on the home page, the user can see the various sections (i.e., workflows manager, jobs manager) appearing on the top. Clicking on them can give access to the other functionalities as discussed below

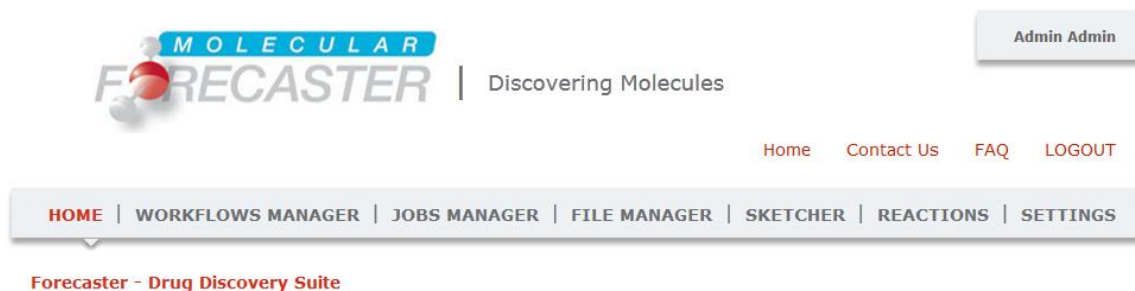

### WORKFLOW MANAGER

The workflow manager shows existing workflows.

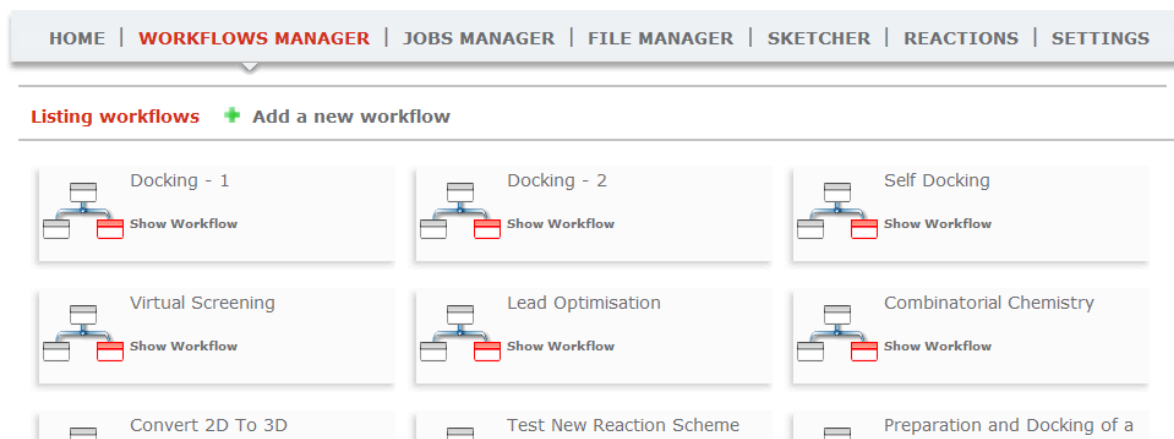

It allows users to view details of existing workflows and to create new ones. Users can only edit and delete workflows that have been created in addition to the out-of-the-box workflows. The workflows can be edited and deleted only when they are not used by any job. In the Figure below Docking 1 is hard-coded and is given as an example by the developers and cannot be edited nor deleted. Workflow 1 was created by a user and can be edited/modified and/or deleted.



To view the schematic details of an existing workflow, click on **Show Workflow**.

## ADDING A NEW WORKFLOW

Clicking on the **+ Add a new workflow** icon allows the user to create a new workflow. A window appears where the user enters a name and a description of the new workflow to be created.

**Add Workflow**


Name

Description


After clicking on **Save**, the user may create a diagram that represents the flow of actions to be executed by clicking on **Start the diagram**

## BUILDING A WORKFLOW DIAGRAM

A new workflow diagram is opened by clicking on **Start the diagram**.

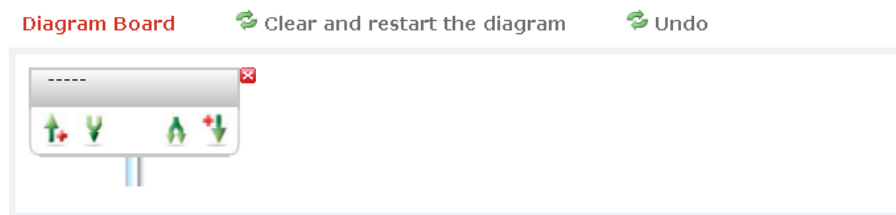


User can add boxes and connectors to the diagram by clicking on the green arrows. A creates a split, creates a merge, and create boxes above and below the current one. By clicking on the icon the user can delete a box from the diagram.

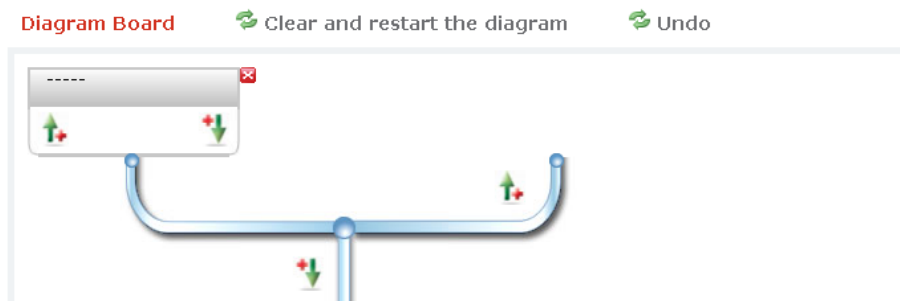

A workflow represents a sequence of actions to be executed. **This sequence runs from left to right and from top to bottom.** In order to create a workflow, the user assigns actions to each




boxes by clicking on the box to be edited and clicking on the  next to one of the actions located on the right hand-side or directly on the action label (i.e., **Add descriptors** ).




Once the diagram is complete, user makes it active and available in the job manager section. The user created workflows become available to all the users on the platform. There is no validation for the compatibility of the actions within the workflow. The user needs to know which actions can be connected together (see the actions section below).

## JOBS MANAGER


The job manager allows users to create a job based on a workflow. The first step is to parameterize the job, the second step is to execute the job, and the final step is to visualize the results.

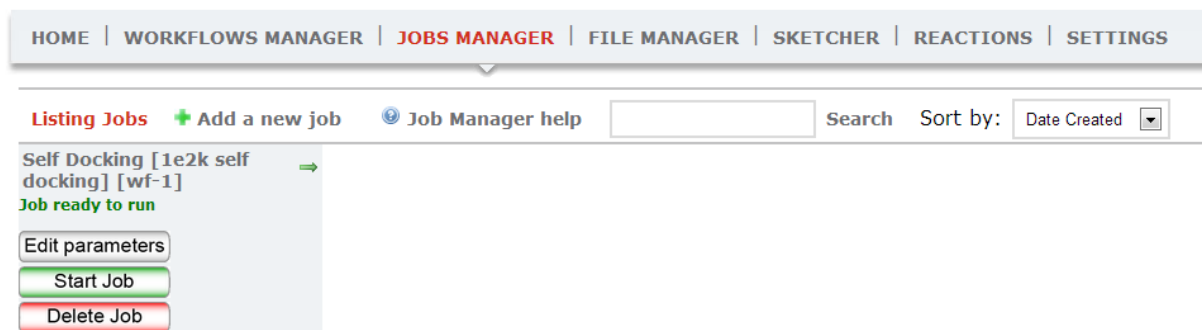



## ADDING A NEW JOB

Users add new jobs by clicking on  **Add a new job** icon, entering the job name, choosing the workflow to be used and clicking continue.

**Add Job**

 Name

 Workflow

 Description

## EDITING THE PARAMETERS AND SELECTING THE FILES TO PROCESS

The editing screen displays the diagram created previously in the workflow manager. In this section, the user has to setup each action of the workflow. This setup is done using a form, which allows users to modify settings and select files to be processed in the job.



**Editing Jobs Parameters**


Name


Workflow

Description


**Current job input files**

  No file chosen 



 [COPY from My Files](#)

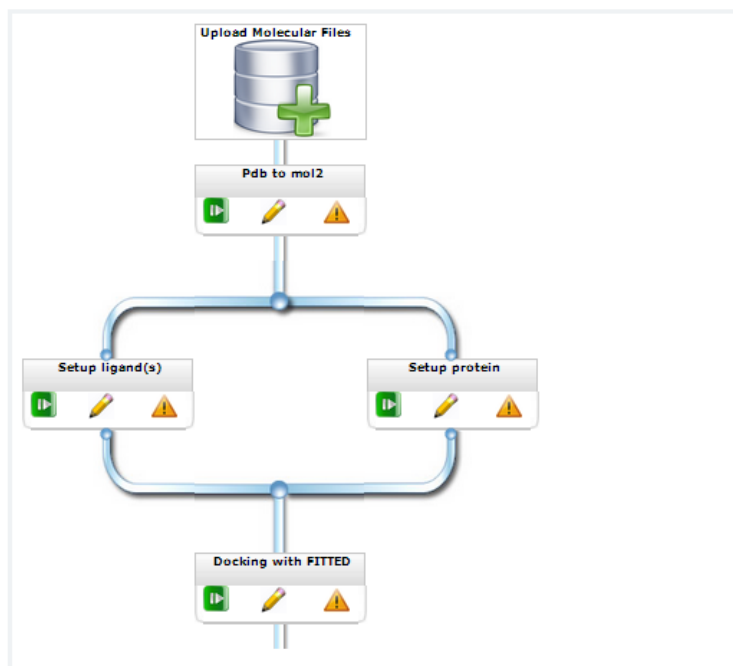**Diagram Board**


First, users have to upload files to be processed by directly uploading them to the system or by copying them from My Files folder. There are two possibilities to add files to the job. Files can be uploaded or transferred from My Files folder using the right-hand side bar.

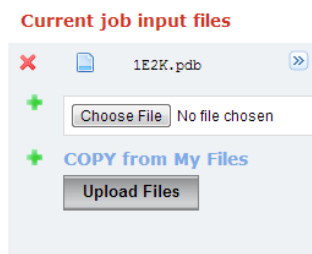


Alternatively, clicking on the container labeled **Upload Molecular Files** brings a window with different possibilities such as drawing directly the molecule, downloading from the web, retrieving from the PDB database or performing a local upload.

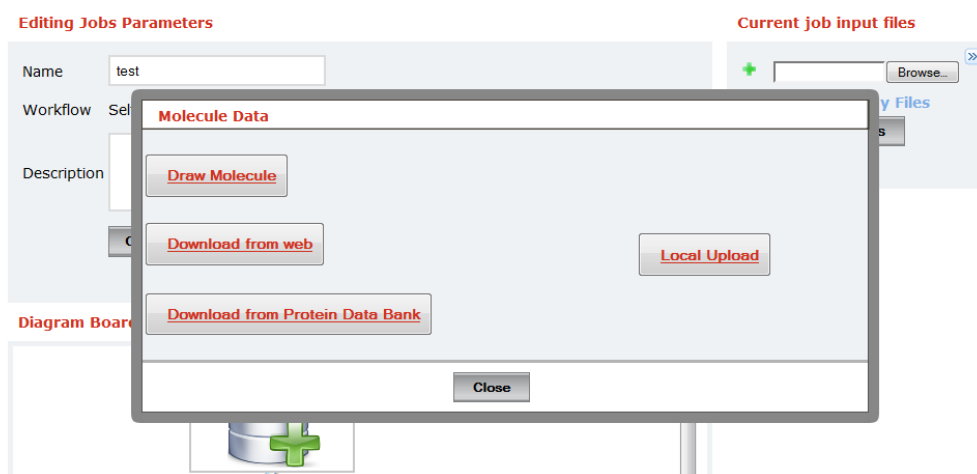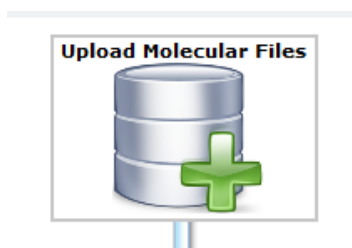

Once the files are uploaded, users can then configure actions. In a workflow, each action (i.e., box) outputs a result file which then becomes the input file for the following action. As the file does not exist yet (the job has not been run yet), an internal predictor guesses the name of these files. As a consequence, **the boxes should be configured from top to bottom.**

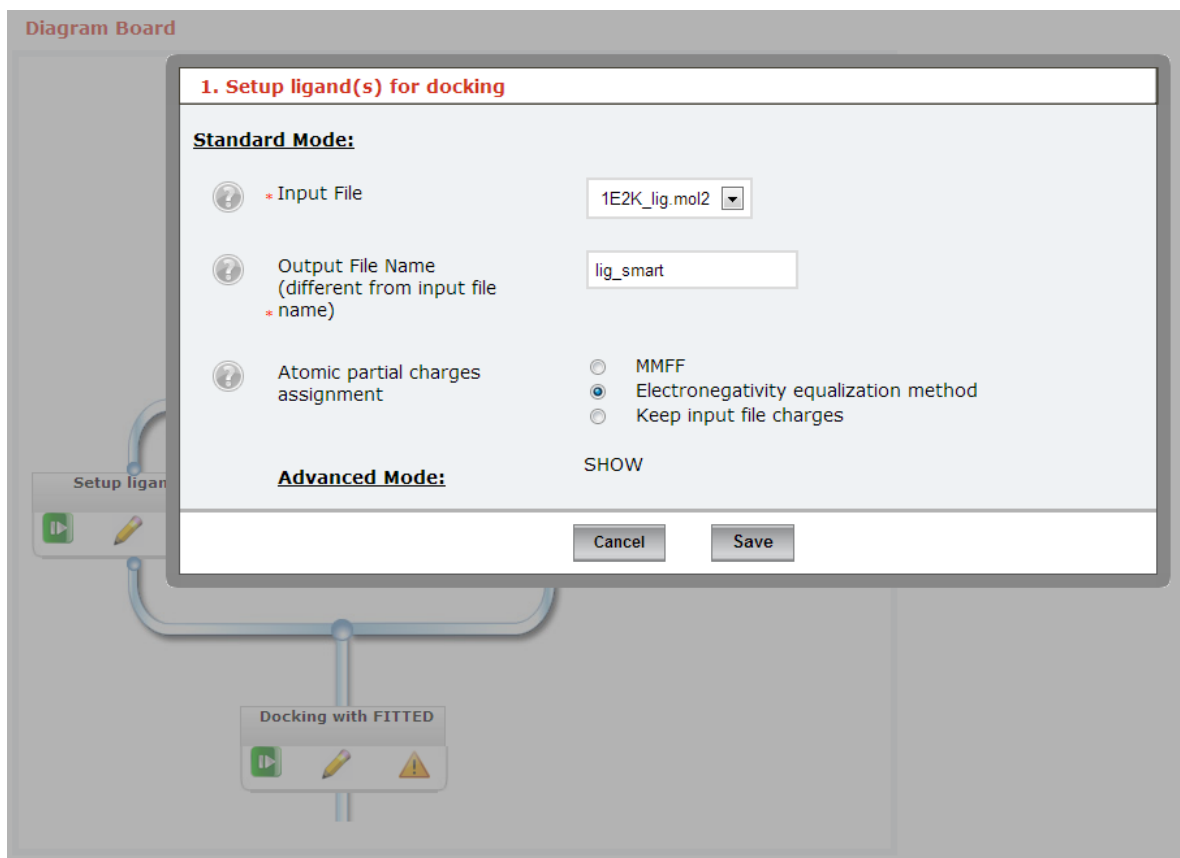

A **pause** feature is implemented which allows the user to stop the execution of a job before a specified action during the workflow so the output of an action can be visualized or edited. The job can then be **resumed** in order to proceed to completion. **This sequence runs from left to right and from top to bottom.** The figure below shows an example where the job will stop after the Prepare protein - pdb to mol2 action from which the protein in mol2 and the ligand in mol2 are generated.

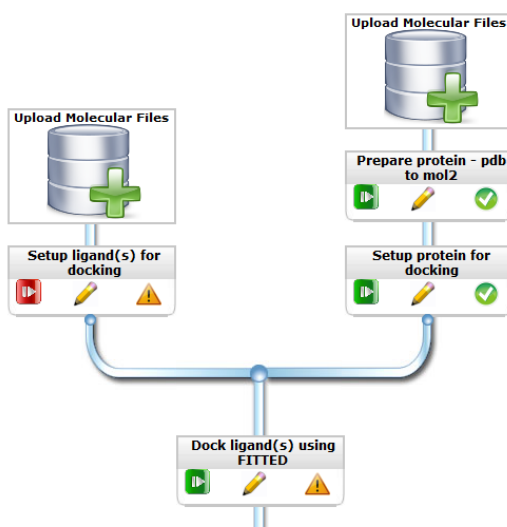

## EXECUTING AND VIEWING JOBS

Once a job is correctly setup, it can be started. The job manager allows users to monitor the job and as each action completes its execution, color led will indicate its status. Output files can be viewed as soon as they are available.

The screenshot shows the Job Manager interface for a job titled "Self Docking [1e2k self docking] [wf-1]". The job status is "Job ready to run". On the left, there are buttons for "Edit parameters", "Start Job", and "Delete Job". The main area lists four steps, each with a red circle icon and a "SHOW PROGRESS" link:

- Prepare protein - pdb to mol2
- Setup ligand(s) for docking
- Setup protein for docking
- Dock ligand(s) using FITTED

On the right, a detailed status panel shows:

- Status : Stand By
- Output : None Available
- Structure : None Available
- Start Date : N/A
- End Date : N/A
- Run Time : N/A

A "close" button is located in the top right corner of the status panel.

## FILE PREVIEW

Once a job is completed (or running), the resulting files can be viewed directly in the job manager, unless the output is too large. In this latter case, the file should be viewed or downloaded in the file manager. A job can be stopped or reset at any time to allow users to re-parameterize and restart it at their convenience. When a job is running, a corresponding folder is created in the file manager where all input and output files are stored. Each job belongs to a particular user and cannot be viewed by other users.

The screenshot shows the Job Manager interface for a job titled "Self Docking [1e2k self docking] [wf-1]". The job status is "Job completed" with a timestamp of "27-06-2013 00:49 AM". There is a "Reset Job" button. The main area lists four steps, each with a green circle icon and a "SHOW PROGRESS" link:

- Prepare protein - pdb to mol2
- Setup ligand(s) for docking
- Setup protein for docking
- Dock ligand(s) using FITTED

On the right, a detailed status panel shows:

- Status : Completed
- Output : docking-results.txt
- Structure : docking\_Complex.pdb
- Start Date : 2013-06-26 at 20:55
- End Date : 2013-06-26 at 20:58
- Run Time : 0h 3min 38sec

A "close" button is located in the top right corner of the status panel. At the top of the interface, there is a "Listing Jobs" section with a search bar and a "Sort by:" dropdown menu set to "Date Created".

## FILE MANAGER

HOME | WORKFLOWS MANAGER | JOBS MANAGER | **FILE MANAGER** | SKETCHER | REACTIONS | SETTINGS

User Files: Upload files to My Files folder File manager help  search

**MY FILES**  
Personal folder  
13 Files

**MY\_PROTEINS**  
5 Files

**TEMP**  
Temporary System Folders  
Empty

**Self Docking [WF-1]**  
1e2k self docking  
Jun 27 2013 00:58:53 10 Files

**Self Docking [WF-2]**  
Self Docking  
Jun 27 2013 00:38:14 1 File

The files in the file manager can be deleted or downloaded. Users can also zip the entire folder or just several selected files. In order to be visualized by the Open Astex Viewer plugin, multiple mol2 files can be combined by clicking on the **COMBINE .mol2** icon. This function assembles several mol2 files (i.e., a protein and a docked ligand) into a single file compatible with Open Astex Viewer. The Icon and List View buttons allow users to change how they view their files.

### FOLDER - MY FILES

Each user has a personal folder called “My Files” in which files from different jobs can be saved to be reused in other jobs. It is used as a transfer folder between a given user’s jobs.

**MY FILES**  
Personal folder  
33 Files

User Files: Upload files to My Files folder File manager help

**SELECT ALL** **ZIP Selected** **COMBINE .mol** **Delete Selected** **COPY to My Files**  
**UNSELECT ALL** **ZIP Folder** **COMBINE .mol2** **Delete Current Folder**

Sort by:  | Filter by:  | Files per page:  | Page:


**MY FILES:**

**2xdw\_lig.mol2**  
[6.2 kb]

**molecule.mol**  
[1.2 kb]

### FOLDER - MY PROTEINS

Each user has a personal folder called “My Proteins” in which proteins files (prepared for docking) from different jobs can be saved to be reused in other docking jobs. It is used inside the “Docking with previously prepared protein files” workflow.

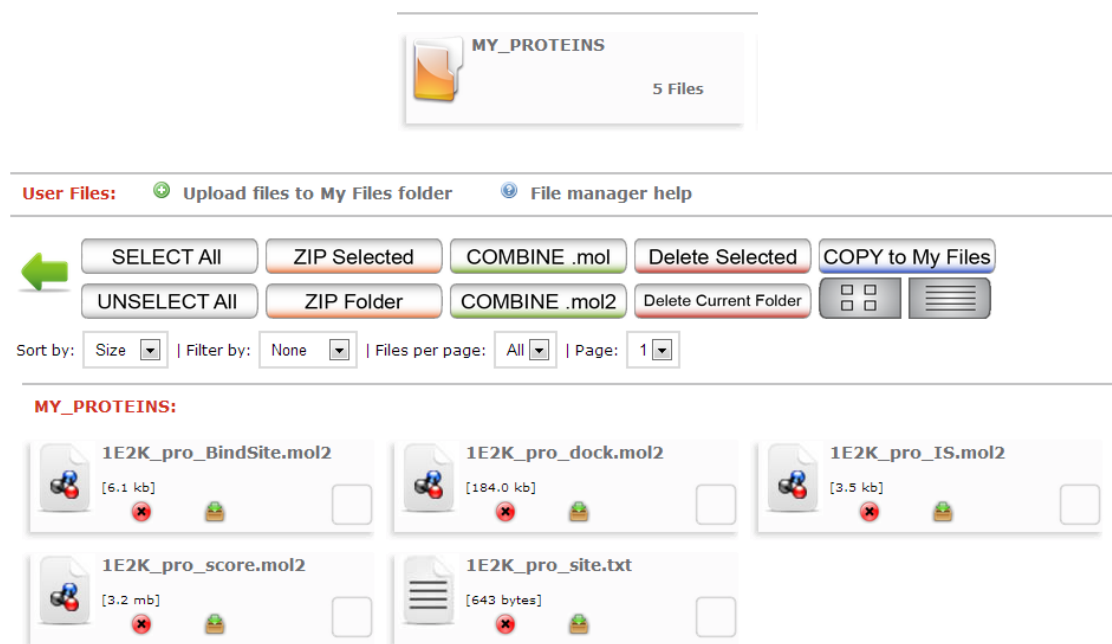

## JOB-SPECIFIC FOLDERS


Files used by each job are saved in the job’s folder using a naming convention such as WF-{id number} e.g. WF-8. The job-specific folder will then read “Not synchronized with the system” and it is at user’s discretion to delete files contained in that folder.



## FILES AND ACTIONS


Files are viewed with different programs depending on their type. All text files are opened in a text editor directly integrated into the interface, while mol2 and pdb files are visualized using the Open Astex Viewer java applet. Files that are bigger than one megabyte (1 Mb) must be downloaded in order to be viewed or edited.

Files located in all folders other than “My Files” folder can be copied to “My Files” by clicking on the “Copy to My Files” icon.

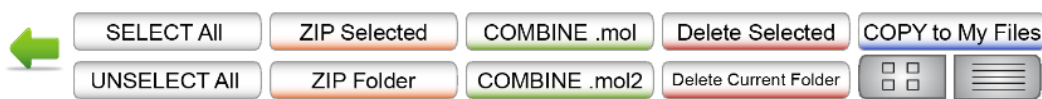


## SKETCHER

The sketcher is a tool that allows users to draw molecules using the SketchEl java applet.

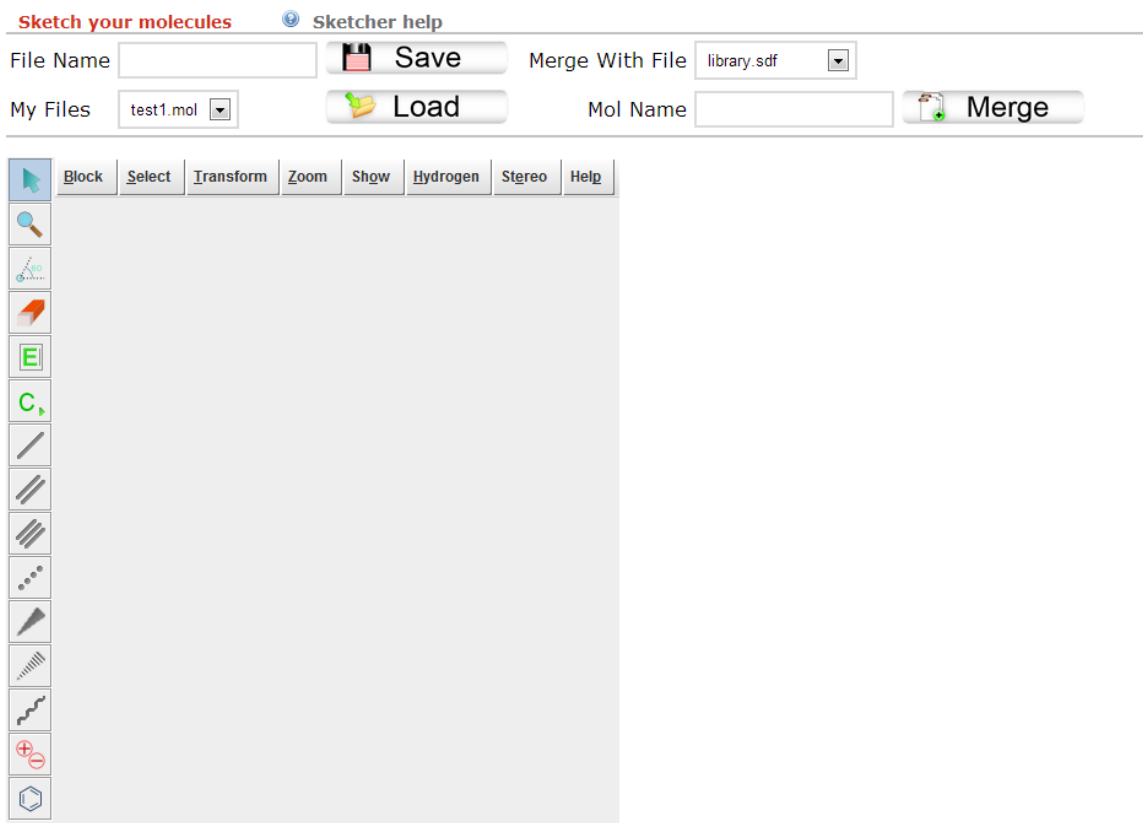

## SAVING AND LOADING FILES

Users can save all molecules they draw using the sketcher in their “My Files” folder. The platform provides an action called “Convert 2D to 3D” which converts the two dimensional drawings (“mol” or “sdf” file format) into three dimensional mol2 files that can subsequently be used as inputs in workflows.

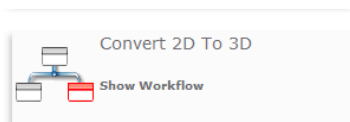

Any molecule that is drawn and saved can later be loaded and reopened in the applet for viewing and editing.




## REACTIONS MANAGER

The reactions manager is a tool that allows users to define reaction rules used in the “Combinatorial Library from 3D Structures” plugin. **This section is in its development phase; we are currently improving the automated reaction rules definitions and developing REACT 2.0 (now in beta version) which will soon fully replace this first version.**

### CREATING REACTIONS


Users have two options when it comes to creating new reaction rules. They can either write the rules themselves or use the SketchEl java applet to draw the reaction.







After assigning a name to the new reaction, user can add as many rules as he likes by clicking on  and filling out the appropriate fields. The rule number defines which reactant the rule applies to, the rule name is usually the reactant functional group, the center is the reacting atom center of the molecule, the kept atoms are the atoms that remain in the product of the reaction, the removed atoms are the atoms that are deleted by the reaction, new atoms are atoms the atoms that are created by the reaction. Once the reaction is defined and saved, it can be used in the Create Combinatorial Library plugin.

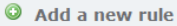HOME | WORKFLOWS MANAGER | JOBS MANAGER | FILE MANAGER | SKETCHER | **SETTINGS**

Editing Reaction

Reaction Name: Amide from carboxylic acid and ar

 Add a new rule


Rule #	Rule Name	Center	Kept Atoms	Removed Atoms	New Atoms	
Rule2	amine	n3	c3, hn	hn		 Remove
Rule1	carboxylic_acid_al	c	c3, o	o		 Remove
Rule1	carboxylic_acid_ar	c	c2, o	o		 Remove
Rule1	carboxylic_acid_vi	c	ca, o	o		 Remove
Rule2	ammonium	n4	c3, hn	hn, hn		 Remove
Rule1	ester	c	c3, o	os		 Remove

Cancel

Save

**\*ADD RULES FROM THE REACTION BELOW**

Alternatively, users can draw and save each reactant and their product using the sketcher, click on Save all to save the drawn structures and have the rules generated automatically by clicking on **\*ADD RULES FROM THE REACTION BELOW**

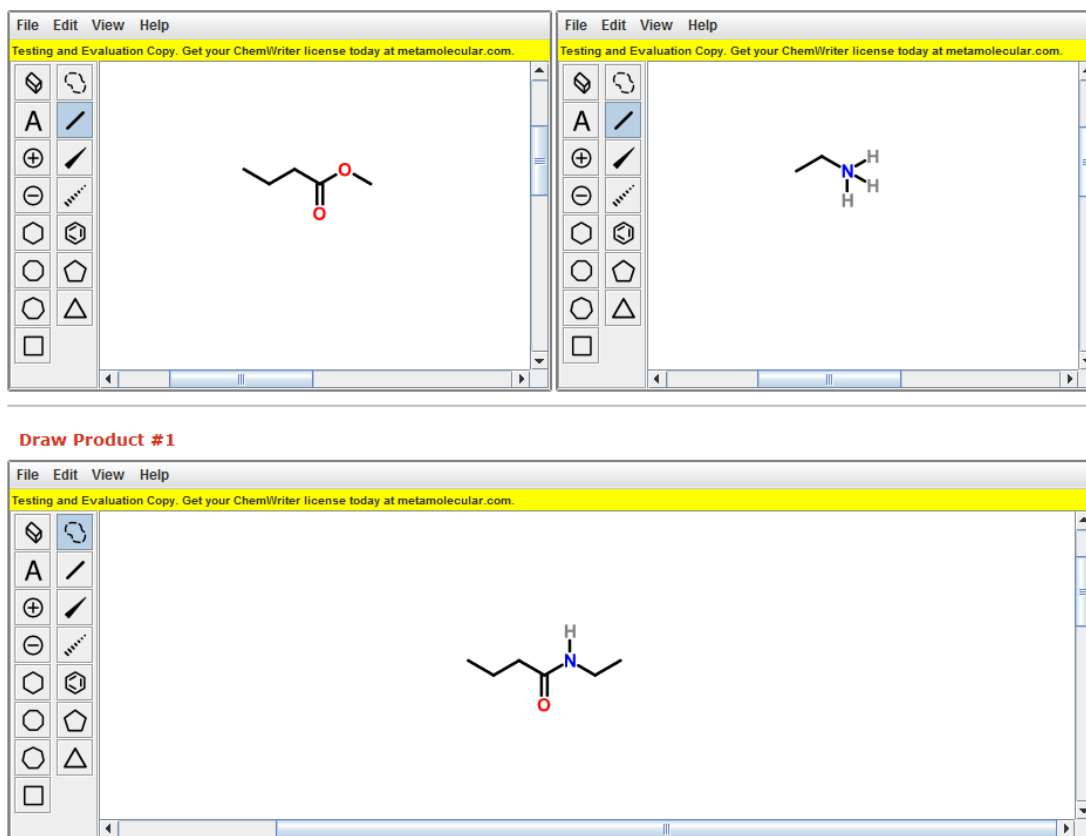


In the current version, the use of the java applet has limitations. Most of the coupling reactions (the most widely used reactions when building libraries) should work. However, rules for transesterification (replacement of an alkoxide group by another alkoxide group) are not figured out by the algorithm currently implemented. Similarly, reactions including two reacting atoms (such as double bond reduction) and metals (Grignard addition) are not yet considered. These issues will be addressed in the next version. If the user cannot produce the rules using the applet, the rules can be set manually. Otherwise, contact Molecular Forecaster ([info@molecularforecaster.com](mailto:info@molecularforecaster.com)) for the development of new rules.


Other limitations are regioselectivity (if more than one functional group of a kind is found, the first one will react) and stereoselectivity (if a chiral center is formed, a single isomer is produced). Stereoselectivity will be considered in the future.


## ADMINISTRATOR INTERFACE - SETTINGS

HOME | WORKFLOWS MANAGER | JOBS MANAGER | FILE MANAGER | SKETCHER | REACTIONS | **SETTINGS**

### Forecaster Advanced Settings

 Plugins Manager

 Users Manager

 Configuration Manager

 Save Database to Fixtures

## USER MANAGER

The user manager is located under the settings menu tab. It allows the administrator to manage users and their roles in the platform.


### USER ROLES AND ACCESS


Each user has a role that is assigned by the platform administrator. Depending on their roles, users are given authorization to access certain features of the program.

### ADDING NEW USERS

Only the program administrator can add or delete users and edit access settings for existing ones.

**Listing Users**  Add a new user

 <div>Admin Admin admin <b>EDIT DELETE</b></div>	 <div>Guest User guest <b>EDIT DELETE</b></div>
---	--

Clicking on the  next to “Add a new user” allows the admin to create a new user. By checking the “Enabled” box, the admin allows the new user to access the application.










### Add User

First name	<input type="text"/>
Last name	<input type="text"/>
Gender	Male <input type="button" value="v"/>
Login	<input type="text"/>
Email	<input type="text"/>
Password	<input type="password"/>
Password confirmation	<input type="password"/>
Enabled	<input checked="" type="checkbox"/>
Admin	<input type="checkbox"/>
<input type="button" value="Cancel"/> <input type="button" value="Save"/>	

## PLUGINS MANAGER

The plugins interface allows users to modify simple information about existing plugins or to delete them.


Listing Plugins   [+ Add a new plugin](#)   [+ New plugin from XML](#)


 Ace 2.0 <a href="#">EDIT</a> <a href="#">DELETE</a>	 Add descriptors 3.0 <a href="#">EDIT</a> <a href="#">DELETE</a>	 Clean structure geometry 3.0 <a href="#">EDIT</a> <a href="#">DELETE</a>
 Convert 2D to 3D 3.0 <a href="#">EDIT</a> <a href="#">DELETE</a>	 Convert file formats 1.0 <a href="#">EDIT</a> <a href="#">DELETE</a>	 Create combinatorial library 1.3 <a href="#">EDIT</a> <a href="#">DELETE</a>
 Dock ligand(s) using FITTED 3.0	 Extract representative library 1.2	 Filter by descriptors 1.3

New actions (programs) can also be integrated into the application. When creating a new plugin, users define all the necessary parameters required for the plugin to work inside of a workflow, such as name, version, and executable files.

## ADDING NEW PLUGINS

Users can add new plugins (i.e., new actions) by clicking on [+ Add a new plugin](#) and defining fields in the plugin interface. In the example described below an action called “convert file format” will be created. This action can next be used in any workflow. This action uses babel as an executable and we would like to create a form to set it up when preparing jobs. The created form is shown below while the preparation of this form and its implementation as an action is given in the subsequent section.

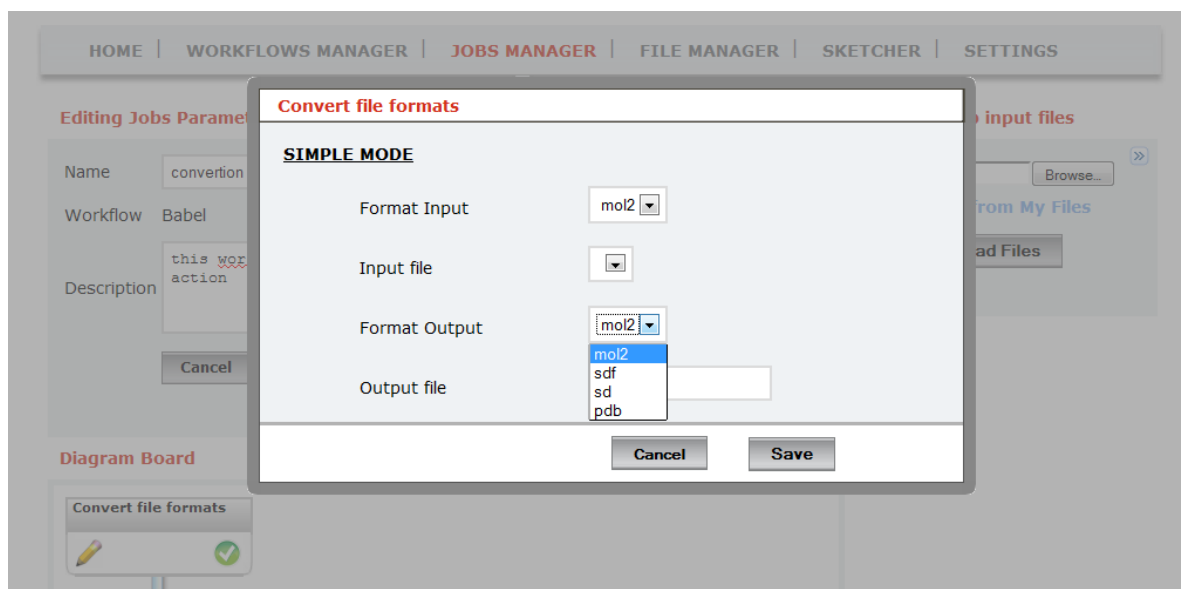

This form is associated to an action (with a name), an executable (i.e., program) and keywords or arguments (if the executable uses a command line). By clicking on  **Add a new plugin** and defining the following window appears and needs to be filled.

Edit Plugin

\* Plugin Name

babel

Unique identifier in the system with no spaces

\* Version

1.0

Release version of the plugin eg. 1.0

\* Function of the plugin

Convert file formats

What the plugin does

Diagram box name

babel

Short description of the plugin function

Active

☒

Displays the plugin on the workflow

Promote

☐

Displays the plugin on the home page

Description

Convert different file formats

Description displayed on the home page

Use system install

☒ Yes

☐ No

Plugin executable path

c:\babel\babel.exe

Use keyword file

☐ Yes

☒ No


Each new plugin has to be assigned a unique name and version. The function of the program may also be described. Additional details about the plugin's purpose may be added in the "diagram box name" field. Users can choose to promote the plugin by checking off the "promote" box, which will ensure that the new plugin is shown on the home page. A description of the plugin's function that will appear on the home page can also be added.

If the plugin's executable file is located on the platform's server, users can use the system install, by checking off "yes" next to "Use system install" and providing a path to the file. If the plugin's executable file is not located on the platform's server, users have to provide it and upload it to the platform.

Some programs work by reading the parameters from a keyword file (usually a text file) while others use various arguments in command line. To instruct the platform to either write a command line or a keyword file, users choose whether or not the arguments for the executable file are located in a keyword file. Selecting "yes" means that arguments are written inside of the keyword file. Selecting "no" means that arguments are specified directly on the command line.

## ADDING PARAMETERS

Users add plugin parameters by clicking on  **Add a new parameter** and filling out the appropriate fields. Each parameter represents a field in the job manager edit screen.

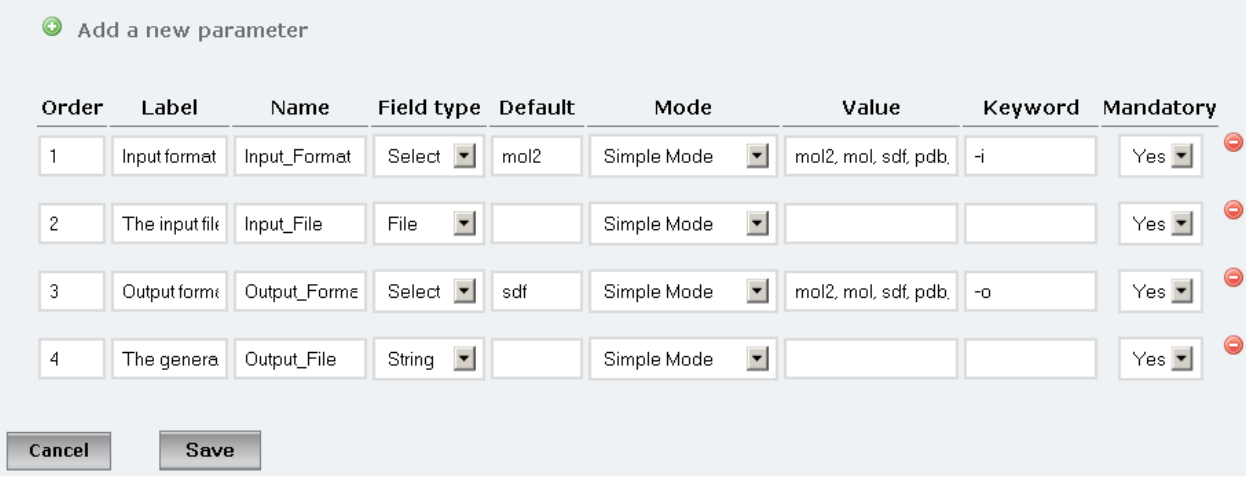

Order	Label	Name	Field type	Default	Mode	Value	Keyword	Mandatory
1	Input format	Input_Format	Select	mol2	Simple Mode	mol2, mol, sdf, pdb,	-i	Yes
2	The input file	Input_File	File		Simple Mode			Yes
3	Output format	Output_Format	Select	sdf	Simple Mode	mol2, mol, sdf, pdb,	-o	Yes
4	The generated file	Output_File	String		Simple Mode			Yes

Cancel Save

**Order** – specify the order in which fields appear on a form

**Label** – describe the field. This label will be used as label in the form

**Name** – specify internal field name in a single word, no spaces.

**Field type** – there are seven possible field types that can be defined.

Check box, radio button or select – when defining one of these field types, the options must be specified in the value field, separated by commas, and the default option will be selected.

String or text – is a text input to be filled out by the user.

File – when this field type is defined, a file is selected from a drop-down menu containing a list of files uploaded to the job using the new plugin.

Hidden – a hidden field contains the keyword value that will be used and that cannot be modified by the user.

Mode – fields can appear in the simple or the advanced mode. The simple mode contains common parameters, modified by most users. Parameters that are less likely to be changed or that are used by more advanced users are placed in the advanced mode.

Keyword – is the name that will be used in the keyword file. It will be followed by the value selected or written in the field, depending on the field type.

Mandatory – specify whether or not filling out a field is mandatory.

Once the form is saved it can be used to edit parameters used in a job.

The screenshot shows a dialog box titled "Convert molecules to different formats" with a sub-header "SIMPLE MODE". It contains four fields, each with a red asterisk indicating it is mandatory:

- Input format of the molecule: dropdown menu with "mol2" selected.
- The input file: dropdown menu with "1hvk\_lig\_1.mol2" selected.
- Output format of the molecule: dropdown menu with "mol2" selected.
- The generated output file: text input field containing "output.mol2".

At the bottom of the dialog are two buttons: "Cancel" and "Save".

## ADDING NEW PLUGINS FROM XML

Users can also create new plugins by uploading and editing their own xml files in the system. They may do so by clicking on **+ New plugin from XML**. Adding a plugin from xml means that users can import an existing or sample definition of a plugin and modify it according to their needs.

The screenshot shows a dialog box titled "Add Plugin from XML". It contains a text input field labeled "XML definition file" with a "Browse..." button to its right. At the bottom of the dialog are two buttons: "Cancel" and "Load XML".


The platform provides a sample file to help in the integration process. The file is located in the doc folder of the application. Once an xml file is loaded, the user gets the new plugin screen, identical to the one that appears when a plugin is created from scratch.


## CONFIGURATIONS MANAGER

The configurations manager can be found in the setting menu tab. It allows the administrator to edit and create new configurations used in the platform.

### DEFAULT SETTINGS


The current settings define the file types that are recognized and can be opened by the application's text editor. The default settings also identify the maximum file size that can be opened by the application without having to be downloaded.

The platform's administrator can add settings by clicking on  **Add a new configuration**.

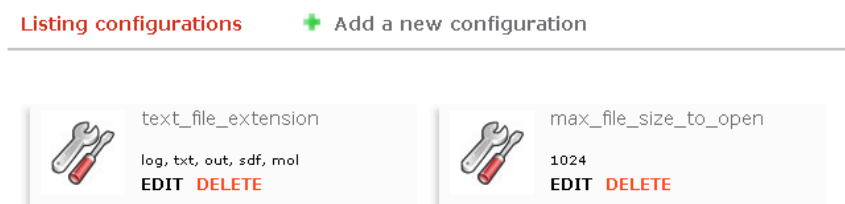





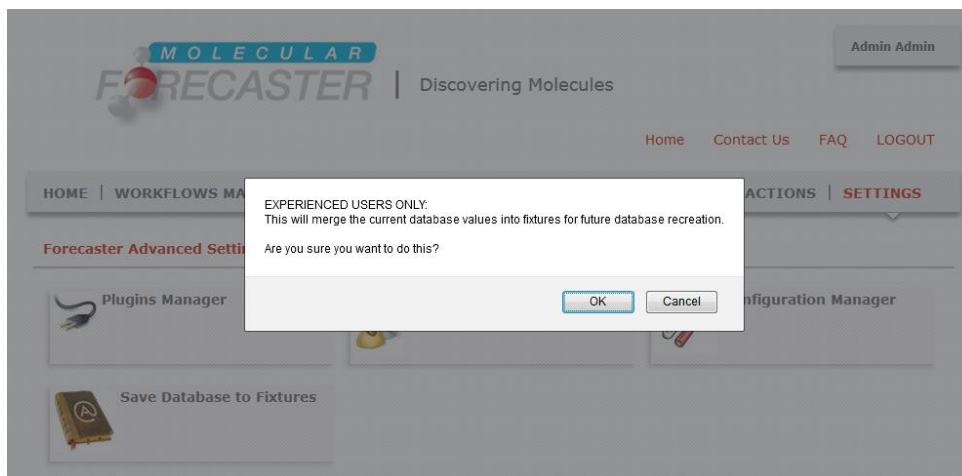

















### SAVE DATABASE TO FIXTURES

This option allows the administrator to store specified database information into files (called fixtures) which is used to re-create the database after an update of the platform. This action saves information about the workflows, jobs, users, plugins, reactions, and configurations that were modified since the first deployment.






List of actions		
<b>Structure-based</b>		
+ 1. dock ligand(s) using fitted		
+ 2. predict som using impacts		
<b>Proteins/Nucleic acids</b>		
+ 1. make proteins similar		
+ 2. superpose protein structures		
+ 3. prepare protein - pdb to mol2		
+ 4. setup protein for docking		
<b>Ligand-based</b>		
+ 1. add descriptors		
+ 2. filter by descriptors		
+ 3. combinatorial libraries from 2d schemes		
+ 4. search for analogues		
+ 5. extract representative library		
+ 6. combinatorial library from 3d structures		
+ 7. add fragments to molecules		
+ 8. pharmacophore identification and search		
+ 9. substructure search		
<b>Small Molecules</b>		
+ 1. setup ligand(s) for docking		
+ 2. add hydrogen to 2d structures		2D
+ 3. convert 2d to 3d		3D
+ 4. clean structure geometry		3D
+ 5. setup reactants for combichem in 3d		
+ 6. structure optimization		
+ 7. add fragments to molecules		2D
<b>Other</b>		
+ 1. convert file formats		
+ 2. function linker		
+ 3. merge files		
+ 4. prepare gamess files		

A number of actions are already integrated in the FORECASTER platform. Details about these actions are given below:

## STRUCTURE-BASED

### DOCK LIGAND(S) USING FITTED

The action Dock ligand(s) using Fitted takes a ligand or a library of ligands and docks it to a protein. It requires that the ligand(s) are setup properly using the action Setup ligand(s) for docking and that the protein has been prepared as well using the actions Prepare protein - pdb to mol2 and Setup protein for docking as shown below. The format of the input files is therefore taken care of by the platform.

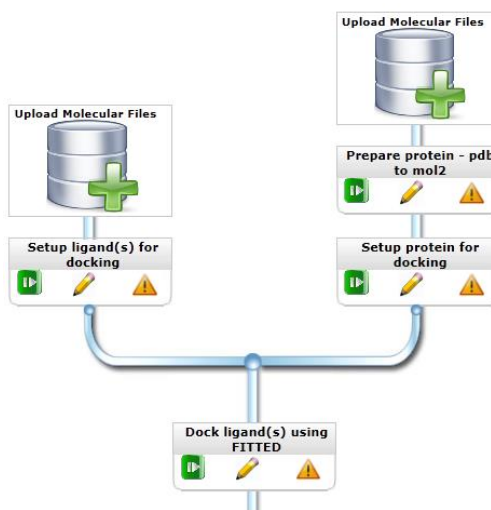

When setting this action in Job manager, a few boxes have to be filled. The number of proteins to be used in the flexible mode is defined in the first box. The protein file name is guessed by the interface as Setup protein for docking should be used in a preceding box.

**Setup protein for docking**

**Standard Mode:**

?
 • Number of protein(s)

1 Protein structure (docking to rigid protein)

?
 Protein File #1

1e2k\_pro.mol2

?
 • Number of Ligand(s)

1 Ligand

?
 Ligand File #1

1e2k\_lig.mol2

?
 Ligand Cutoff

7

?
 Output File Name  
(different from input file name)
 • name)

process

?
 • Prepare for

Docking to flexible protein

**Advanced Mode:**

SHOW

Cancel Save

## PREDICT SOM USING IMPACTS

This action uses the IMPACTS program and predicts the most likely site(s) of metabolism (SoM) and transition state (TS) structures of small molecules when reacting with the CYP heme as the activated iron-oxygen species. It uses a significantly modified version of our FITTED docking program to predict CYP-mediated metabolism of small molecules.

The action contains a single box, the molecule can be drawn in 2D within the sketcher (or uploaded as an sdf file). A 3D format is also supported (single or multi-mol2 files). The user can choose between the five implemented CYP 450 namely: 1A2, 2D6, 2C9, 2C19 and 3A4 or all five in a single run.

The screenshot shows a software window titled "2. Predict SOM using IMPACTS". It has a "Standard Mode:" section and an "Input / Output parameters" area. There are three rows of parameters, each with a question mark icon on the left. The first row is "P450 enzyme" with a dropdown menu showing "2C9". The second row is "Ligand File" with a dropdown menu showing "flurbiprofen.mol". The third row is "Output File Name (different from input file name)" with a text box containing "impacts\_results". Below these parameters is a "SHOW" button. At the bottom of the window are "Cancel" and "Save" buttons. There is also an "Advanced Mode:" section which is currently collapsed.

## PROTEINS/NUCLEIC ACIDS

### MAKE PROTEINS SIMILAR

This action takes multiple protein pdb files (the proteins need to have greater than ~90% sequence identity) and perform a sequence alignment to superpose all the proteins together. In a second step, all the proteins are made similar by deleting or mutating amino acids to be identical to the first in the list (the reference protein). Corresponding pdb files are then obtained and can be used with the action *Prepare pdb to mol2*. If the sequence identify is too low, an error message will be given and no structure produced.

**1. Make proteins similar**

**Standard Mode:**

Number of protein(s) 3 Proteins

Protein File #1 1E2K.pdb

Chain(s) All

Protein File #2 1KIM.pdb

Chain(s) All

Protein File #3 1KI3.pdb

Chain(s) All

Output

Identify ligand residues in pdb 1KI3.pdb

Number of ligand residues 3

	Residue Name	Chain Name	Residue Number
Residue 1	TMC	A	500
Residue 2	THM	A	1
Residue 3	PE2	A	1

Cancel Save

## SUPERPOSE PROTEIN STRUCTURES

This action takes multiple protein structures (pdb files) and performs a sequence alignment in order to superpose all the proteins together. Corresponding pdb files are then obtained. This action is used in the “Cross docking” workflow.

**2. Superpose protein structures**

**Standard Mode:**

Number of protein(s) 2 Proteins

Protein File #1 1E2K.pdb

Chain(s) All

Protein File #2 1KIM.pdb

Chain(s) All

Output superpose\_output

Superposition all pairs No

Number of ligand residues 2

	Residue Name	Chain Name	Residue Number
Residue 1	TMC	A	500
Residue 2	THM	A	1

Cancel Save

## PREPARE PROTEIN – PDB TO MOL2

This action takes a protein pdb file, add hydrogens to the protein (according to the residue pKa), and to water oxygen atoms, searches for the optimal rotamers for asparagines, glutamines and histidines, reconstructs and optimizes missing side chains, extract the ligand and finally outputs a protein file and a ligand file in mol2 format. This format is appropriate for most programs and required by most actions in FORECASTER as shown below.

**3. Prepare protein - pdb to mol2**

Standard Mode:

Number of protein(s) 1 Protein

Protein File #1 1E2K.pdb

Output prepare\_protein

Identify ligand residues in pdb Select pdb

Number of ligand residues 1

	Residue Name	Chain Name	Residue Number
Residue 1	TMC	A	500

Re-assign hybridization No

Protonate atom No

Optimize ☒ Yes ☐ No

Iterations 5

Side-chain conformations ☐ Generate new side chain conformations ☒ Take from input file only


Water Molecules Crystallographic

Macromolecule Protein

Advanced Mode: SHOW

Cancel Save

When setting this action in job manager, the user will have to provide the ligand name (i.e., TMC A 500 below). This information can be found in the pdb file. By selecting the desired pdb file in the “Identify ligand residues in pdb” section, Open Astex viewer will open and allow the user to click on the ligand and obtain the ligand residues code.

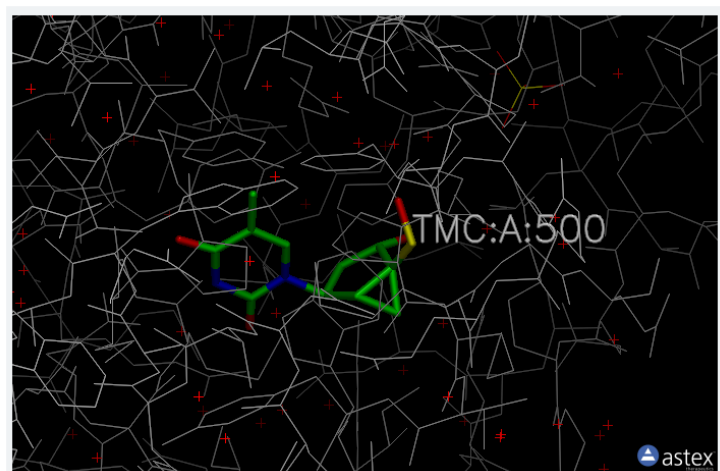

## SETUP PROTEIN FOR DOCKING

This action is required to prepare the necessary files for FTTED to work. In the publications describing FTTED [1-3], we refer to the use of a program PROCESS which is the core of this action. It requires a protein mol2 file as input. Providing a ligand mol2 file helps identifying the binding site. The ligand file can be obtained by the action `Prepare protein - pdb to mol2`.

**4. Setup protein for docking**

**Standard Mode:**

- Number of protein(s): 1 Protein structure (docking to rigid protein)
- Protein File #1: 1E2K\_pro.mol2
- Macromolecule: Protein
- Number of Ligand(s): 1 Ligand
- Ligand File #1: 1E2K\_pro.mol2
- Ligand Cutoff: 7
- Output File Name (different from input file name): process\_protein
- Prepare for: Docking to flexible protein
- Keep files for later use: ☒ Yes ☐ No

**Advanced Mode:** SHOW

Cancel Save

## LIGAND-BASED

### ADD DESCRIPTORS

As described in ref. 2, SMART can add descriptors that can be used for further filtering. This action has been built from this program. The ligand must be in 3D and have hydrogen atoms added when given to this action. Formats such as sdf and mol2 are accepted and automatically detected by the action. Within the `add descriptors` action, a number of descriptors and functional groups are automatically identified and added as a bitstring in the resulting mol2 output file.

**1. Add descriptors**

**Standard Mode:**

Input File:

Input Format: ☒ 3D - sdf or mol2 ☐ 2D - sdf

Output File Name (different from input file name):

**Advanced Mode:**

## FILTER BY DESCRIPTORS

Filter by descriptors uses the descriptors prepared using the action mentioned above and keep only the molecules with the appropriate properties. It uses the REDUCE program.

Descriptor	Minimum	Maximum
Molecular weight	<input type="text" value="0"/>	<input type="text" value="500000"/>
Net total charge	<input type="text" value="-20"/>	<input type="text" value="20"/>
Number of hydrogen bond acceptor(s)	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of hydrogen bond donor(s)	<input type="text" value="0"/>	<input type="text" value="500"/>
Total number of atoms	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of heteroatom(s)	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of oxygen atom(s)	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of nitrogen atom(s)	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of sulphur atom(s)	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of metal atom(s)	<input type="text" value="0"/>	<input type="text" value="0"/>
Rings	<input type="text" value="0"/>	<input type="text" value="1000"/>
Number of rotatable bond(s)	<input type="text" value="0"/>	<input type="text" value="6000"/>
Ionizable group(s)	<input type="text" value="0"/>	<input type="text" value="200"/>
logP	<input type="text" value="0"/>	<input type="text" value="5"/>
Polar Surface Area	<input type="text" value="0"/>	<input type="text" value="500"/>
logS (Solubility)	<input type="text" value="0"/>	<input type="text" value="500"/>


  

Functional group	None	Filter	Optional	Minimum	Maximum
Acyl chloride	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>
Aldehyde	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>
Alkyl chloride	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>
Alkyl bromide	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>
Alkyl iodide	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>
Amide	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="text" value="0"/>	<input type="text" value="0"/>

## COMBINATORIAL LIBRARY FROM 2D STRUCTURES

This action is based on our program REACT 2.0 that takes two libraries of chemicals (in 2D) and a reaction scheme (either from three separate mol files or a MDL rxn scheme file) to prepare a combinatorial library. In this new action, the reaction is not defined using the

reactions manager. The reaction scheme is given to the interface from the dialog box along with the two libraries of chemicals. To define the reaction, the reactant #1, reactant #2 and the product files in MDL mol format are required.

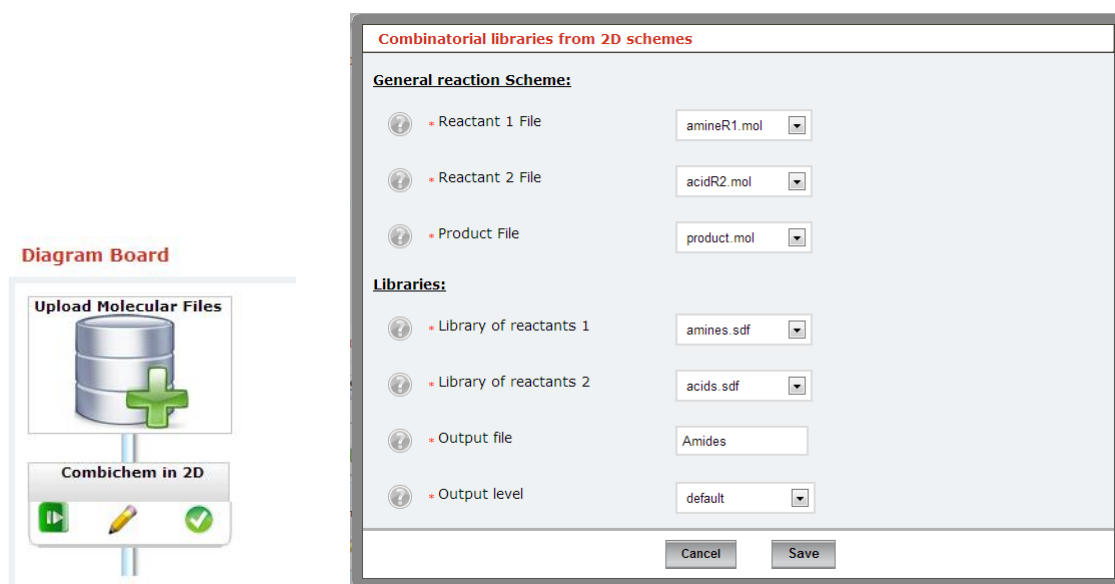

---

## SEARCH FOR ANALOGUES


This action takes a library and a hit molecule and creates a new library of molecules from the library that are similar to the hit. It uses the program SELECT. It can take either 2D or 3D structure files as input.

---

## EXTRACT REPRESENTATIVE LIBRARY

This action takes a library (either 2D or 3D), clusters the molecules by similarity and outputs a smaller library. It uses the program SELECT. The clustering allows the user to remove similar molecules (e.g., to keep the same diversity while reducing the number of molecules). The number of clusters is the maximum number of molecules extracted based on the value of the Tanimoto coefficient used.



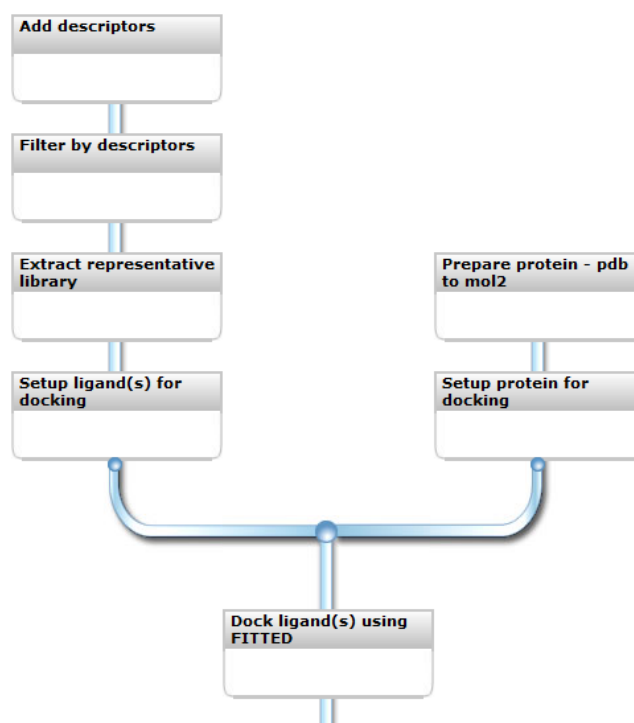

**Extract representative library**

**Standard Mode:**

Input Library File:

Output File Name (different from input file name):

Number of Clusters:

Minimum Tanimoto coefficient (0-100):

Format of input file: ☒ Standard MOL2 ☐ 3D SD or 3D SDF ☐ FORECASTER

**Advanced Mode:**

## COMBINATORIAL LIBRARY FROM 3D STRUCTURES


This action is based on our program REACT 1.0 that takes two libraries of chemicals (in 3D only) and a reaction scheme to prepare a combinatorial library. The reactants name are guessed by the interface as `Setup reactants for combichem` should be used in a preceding box as shown in the Figure below. The reaction is defined and pre-encoded via the reactions manager.

**Create combinatorial library**

**Standard Mode:**

- Input Reactant 1: scaffold-1.mol2
- Input Reactant 2: scaffold-1.mol2
- Reaction type: Suzuki coupling

Cancel Save

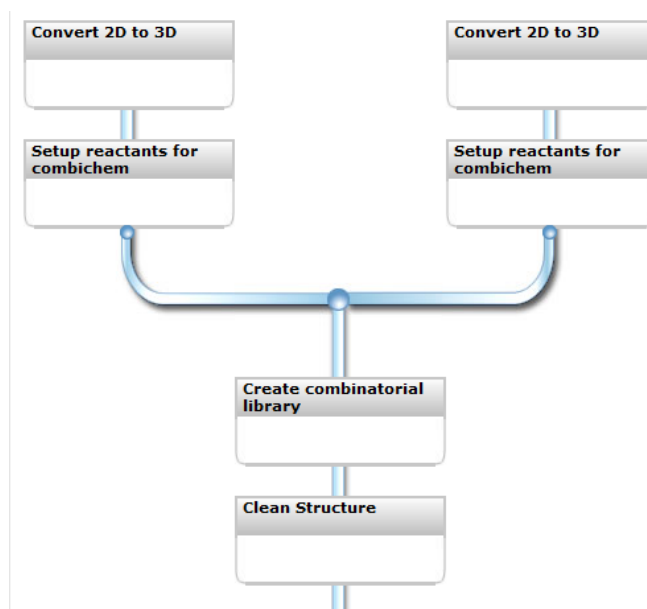

## ADD FRAGMENTS TO MOLECULES

This action is in development and not fully functional yet. Contact us for more information.

## PHARMACOPHORE IDENTIFICATION AND SEARCH

This action uses the RESHAPES program which is a tool for 3D pharmacophore search. It derives a set of weighted pharmacophores from known actives and uses these pharmacophores to identify potential active molecules in a library. Both the actives and the library files need to be in 3D sdf or mol2 formats. **This new program is still in development (beta).**

**8. Pharmacophore identification and search**


**Standard Mode:**

- Actives: reshape\_test\_actives.sdf
- Library to be screened: reshape\_test\_library.sdf
- Output File Name (different from input file name): potential\_actives

Cancel Save

## SUBSTRUCTURE SEARCH


This action uses FINDERS to search for molecules featuring a given substructure within a database of chemicals. The database format is a normal sdf file without any preparation. This action is still in development and not yet fully validated.

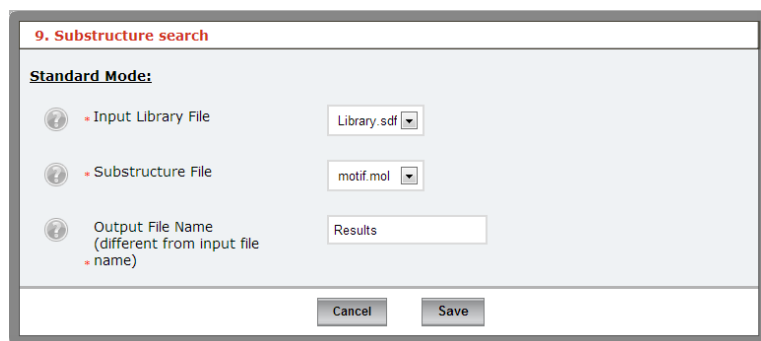

## SMALL MOLECULES

### SETUP LIGAND(S) FOR DOCKING.

In the publications describing FITTED [1-3], we referred to SMART as a program to setup the ligand files. This action has been built from this program. The ligand must have hydrogen atoms added when given to this action (the user may use Convert 2D to 3D to add hydrogens or obtain this ligand from the Prepare protein - pdb to mol2 action). 3D structure files in sdf and mol2 formats are accepted and automatically detected by the action.

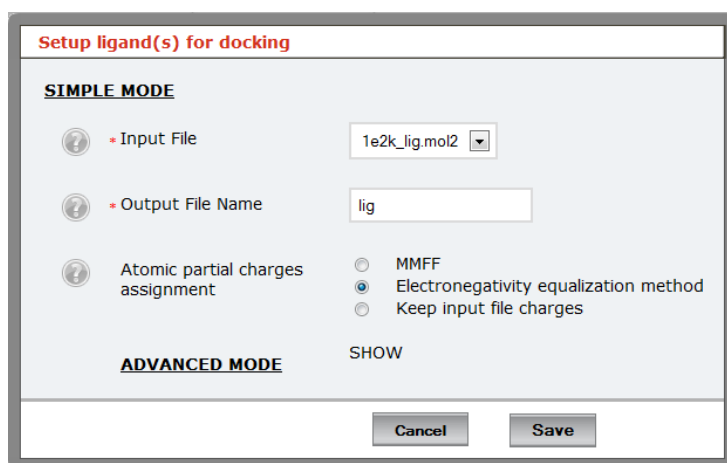


### ADD HYDROGEN TO 2D STRUCTURES

This action adds hydrogens to 2D molecules. It is used by different actions when 3D is not required but hydrogens need to be present.

### CONVERT 2D TO 3D

This action takes the input structure(s), adds missing hydrogens, converts from 2D to 3D and optimizes through energy minimization. The file formats compatible with this action are .mol and .sdf (2D).

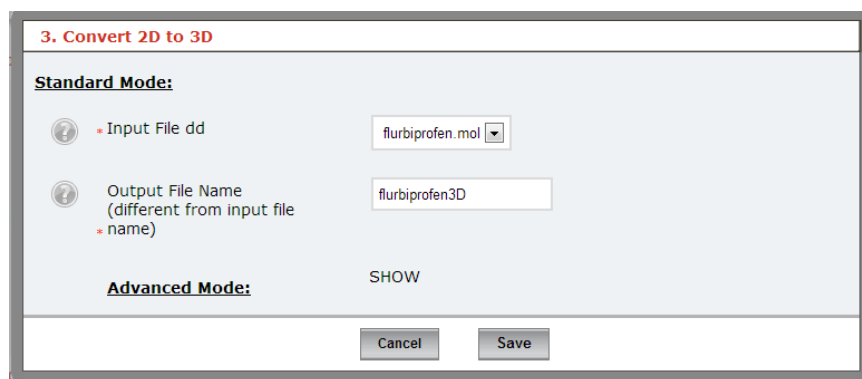

---

### CLEAN STRUCTURE GEOMETRY

This action cleans (add missing hydrogens, ...) and optimizes the input structure(s) through energy minimization (optional). It requires 3D structure files (mol2 or sdf).

---

### SETUP REACTANTS FOR COMBICHEM IN 3D

This action prepares the reactants for further processing with the program REACT implemented in the action below. The dialog box is similar to the action `add descriptors`. It requires 3D structure files (mol2 or sdf).

---

### STRUCTURE OPTIMIZATION

This action optimizes structures that have been setup for docking (conjugate gradient minimization). It requires 3D structure files (mol2 or sdf).

### OTHER

---


### CONVERT FILE FORMATS

This action uses babel to convert a structure from one format (i.e., mol2) to another one (i.e., sdf).

---

### FUNCTION LINKER

This is not a real action. This action adds a connector (i.e., a vertical bar) to the workflow.

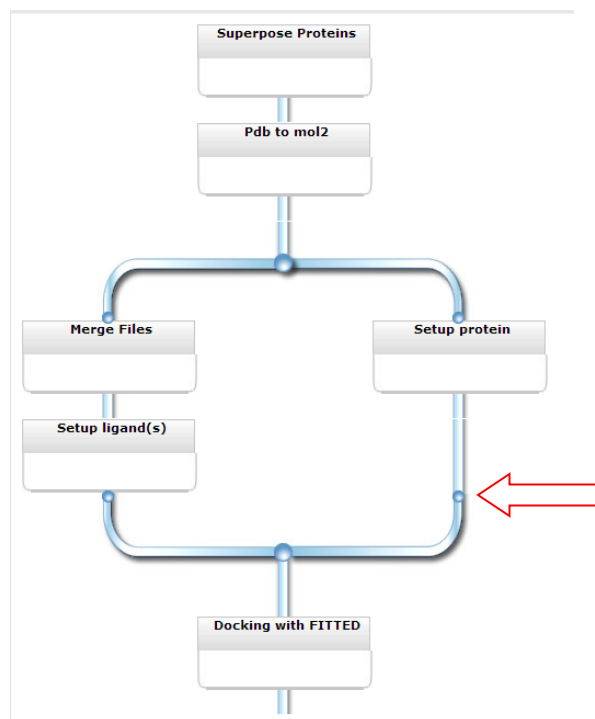

## REFERENCES

1. Christopher R. Corbeil, Pablo Englebienne, and Nicolas Moitessier. Docking Ligands into Flexible and Solvated Macromolecules. 1. Development and Validation of FITTED 1.0. *J. Chem. Inf. Model.*, **2007**, 47 (2), pp 435–449.
2. Christopher R. Corbeil, Pablo Englebienne, Constantin G. Yannopoulos, Laval Chan, Sanjoy K. Das, Darius Bilimoria, Lucille L’Heureux and Nicolas Moitessier. Docking Ligands into Flexible and Solvated Macromolecules. 2. Development and Application of FITTED 1.5 to the Virtual Screening of Potential HCV Polymerase Inhibitors. *J. Chem. Inf. Model.*, **2008**, 48 (4), pp 902–909.
3. Christopher R. Corbeil and Nicolas Moitessier. Docking Ligands into Flexible and Solvated Macromolecules. 3. Impact of Input Ligand Conformation, Protein Flexibility, and Water Molecules on the Accuracy of Docking Programs. *J. Chem. Inf. Model.*, **2009**, 49 (4), pp 997–1009
4. Pablo Englebienne and Nicolas Moitessier. Docking Ligands into Flexible and Solvated Macromolecules. 4. Are Popular Scoring Functions Accurate for this Class of Proteins? *J. Chem. Inf. Model.*, **2009**, 49 (6), pp 1568–1580
5. Roxanne Kieltyka, Pablo Englebienne, Johans Fakhoury, Chantal Autexier, Nicolas Moitessier and Hanadi F. Sleiman. A Platinum Supramolecular Square as an Effective G-Quadruplex Binder and Telomerase Inhibitor. *J. Am. Chem. Soc.*, **2008**, 130 (31), pp 10040–10041

6. Pablo Englebienne and Nicolas Moitessier. Docking Ligands into Flexible and Solvated Macromolecules. 5. Force-Field-Based Prediction of Binding Affinities of Ligands to Proteins. *J. Chem. Inf. Model.*, **2009**, 49, pp 2564–2571
7. Janice Lawandi, Sylvestre Toumieux, Valentine Seyer, Philip Campbell, Sabine Thielges, Lucienne Juillerat-Jeanneret and Nicolas Moitessier. Constrained Peptidomimetics Reveal Detailed Geometric Requirements of Covalent Prolyl Oligopeptidase Inhibitors. *J. Med. Chem.*, **2009**, 52 (21), pp 6672–6684